



**GB      Instruction Manual Item No. LH-538B**  
**Washable Heating Pad**  
**Please read before using!**

## Introduction




Thank you for choosing the Washable Heating Pad. The Heating Pad is a product of high standard that will last a long time if properly used and handled. There are four heat settings level, the Heating Pad has a safety switch that will switch off automatically after 120 minutes.

The heat is regulated to ensure a permanent even heat distribution. If in use you will not experience the changes between on and off, or heating up and cooling down. A sensor cable provides for double protection against overheating.

The product is tested according to European standard and therefore technically safe. Please read the safety advice in your user manual.

## Washable Heating Pad Components

1. Heating Pad
2. High Heat with red LED indicator
3. Medium Heat with orange LED indicator
4. Low Heat with yellow LED indicator
5. Warm Heat with green LED indicator
6. ON/OFF button
7. Connection clutch
8. Controller

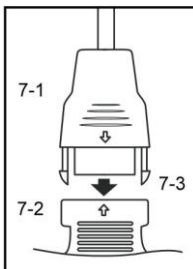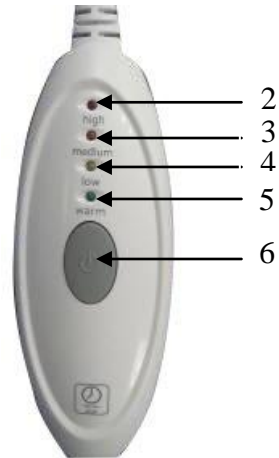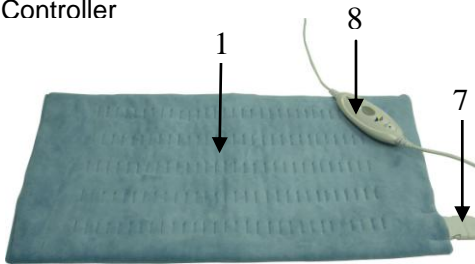

- 7-1. Controller's Connector  
7-2. Wrap's Connector  
7-3. Locking Clips

**The appliance is not to be used by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction. Children being supervised not to play with the appliance.**

## **HOW TO USE**

Heating Pad is specifically designed to fit the waist, back, knee and everywhere you like to use on your body, where users experience the most physical signs of tension and stress. Heating Pad comfortably on area of your body where you like.

Make sure Connection Clutch (7) is securely connected.

The Heating Pad is used for well-being and not for any medical purposes. Heat influences the blood flow and relaxes the muscles. The Heating Pad itself can be removed for washing.

### **Turn on heat**

Put plug into 220V AC electrical outlet. The Heating Pad is ready for use. Select from 4 heat settings by pressing ON/OFF button successively to achieve desired heat level. To turn off, you must scroll through the remaining heat settings by pressing the ON/OFF button sequentially. The 4 heat settings with different LED indicators:

- The green LED indicated for Warm level.
- The yellow LED indicated for Lower level.
- The orange LED indicated for Medium level.
- The red LED indicated for High level.

Note: If the heat level indicators do not light up or stay on, the plug or control connector may not be firmly and completely connected. Unplug from wall outlet, check the Heating Pad connection again and re-plug into outlet.


### **Automatic switch off**

120 minutes after switching on, the Heating Pad will switch off automatically. If the Heating Pad should stay in use, you have to press ON/OFF button again.

**Attention:** Extended use may lead to burned skin.

## Maintenance and Storage

The Heating Pad can be washed in the washing machine after removing the Controller. Remove the mains cable at the connection clutch (7). You may return the Heating Pad either to its shipping carton or store flat in a cool, dry place. Take care not to damage the surface of the Relax Wrap with sharp objects. The control unit must not be allowed to wash.

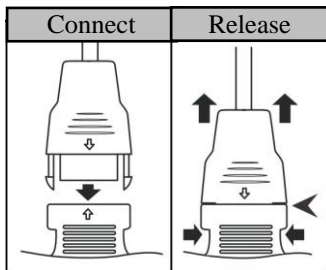

## Symbol explanation



Please, read the instructions



Do not use if folded or creased



Do not insert pins or needles

## Important Safety hints

# RETAIN MANUAL FOR FURTHER USE



- Read the user manual thoroughly
- Retain the manual for future use
- Do not prick in with needles or other sharp objects.
- Heating Pad
  - only plug into adequate power point
  - do not use if wet or damp
  - do not use on heat insensitive persons, children or pets
  - do not use if bent or folded
  - do not use on helpless persons
  - do not fall asleep if Heating Pad is switched on
  - do not fold Heating Pad if still warm
- do not place heavy or edged objects on the Heating Pad, which may damage or bend the electrical wires inside the Heating Pad (eg. Suit cases )
- Keep plug and connecting cable dry at all times.
- The switch will get warm and therefore should not be covered, danger of overheating.
- Persons who only experience heat to a certain level should only use the Heating Pad under supervision
- If the Heating Pad or parts of it are damaged the Heating Pad should not be used
- After use pull plug out off power point
- Do not place other heaters or heat producing objects onto the switched on cushion.
- Always supervise if used for children
- Avoid creasing the appliance
- Examine the appliance frequently for signs of wear or damage. If there are such signs or if the appliance has been misused, return it to the supplier prior to any further use
- This appliance is not intended for medical use in hospitals
- If the supply cords damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- Do not use if wet.

- When storing the appliance, allow it to cool down before folding.
- Do not crease the appliance by placing items on top of it during storage.
- When the heating cushion is likely to be used for an extended period, for instance by the user falling asleep, controls should be adjusted to a setting recommended for continuous use.
- The detachable control unit is only to be used with the types that are marked on the appliance.

## Delivery contents

- Heating Pad
- User manual

## Technical data

Model:	LH-538B
Alternating current	220-240V~ 50Hz
Max. power	ca. 100W
Measurements	300*600mm



**For the sake of the environment**, we request that you take packaging material to be recycled (collection points).



## Disposal of appliances and batteries

After the implementation of the European Directive 2002/96/EU in the national legal system, the following applies:

Electrical and electronic devices may not be disposed of with domestic waste. Consumers are obliged by law to return electrical and electronic devices at the end of their service lives to the public collecting points set up for this purpose or point of sale. Details to this are defined by the national law of the respective country. This symbol on the product, the instruction manual or the package indicates that a product is subject to these regulations. By recycling, reusing the materials or other forms of utilising old devices, you are making an important contribution to protecting our environment.

## **Guarantee**

The guarantee period is 12 months from the date of purchase and covers all manufacturers' errors in material and quality. This guarantee only applies if the instructions for use have been followed and becomes ineffective if force has been applied to the appliance or it has been used in any but the correct and proper manner, or if it has been serviced by any unauthorised person.