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
English

1 Introduction

Thank you for choosing the SGL Rechargeable Warming Belt. This product is designed to be an auxiliary source of warmth for individual use that is unencumbered by wires tied to a power source.

2 The Components

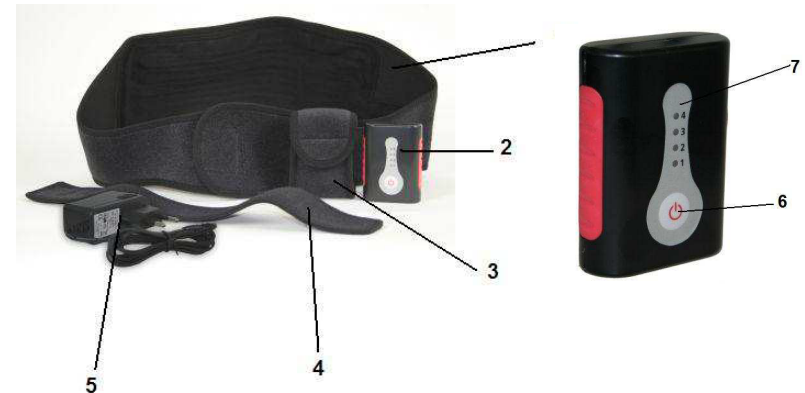- 1) Rechargeable Warming Belt
- 2) Controlling electronics and Lithium Ion rechargeable battery
- 3) Pocket for controller with hook and loop closure
- 4) Belt extender
- 5) Charger
- 6) Start/Stop key
- 7) Light Emitting Diodes (LEDs)

3 Description of functions

The “Rechargeable Warming Belt” is supplied with power by a Lithium Ion rechargeable battery that is integrated in the device’s controller. Simply charge the battery by plugging into power main. After being charged, the battery/controller piece is plugged into the belt and can be carried in the pocket attached to the belt by a hook/loop attachment. The degree of warmth can be selected in 4 levels between 38 and 55 degrees Centigrade.

The battery/controller component of this device can be procured under the part number LH-115.30 in case a second battery is desired.

4 Use

Please charge the battery fully before using the device for the first time. To accomplish this, plug the charger (5) into a power main and connect with the controller/battery (2) and allow charging for 3-4 hours. The extent of the charge is indicated by a green LED on the controller.

Afterwards:

- a) Connect the cable of the warming belt with the controller.
- b) Depress the Start/Stop key for 2 seconds
- c) Select the desired level of warmth by pressing the Start/Stop key repeatedly.
- d) Place the controller (2) in the side pocket during use. When doing this, take care not to allow access to the Start/Stop key during use.
- e) Attach carrying pocket to the hook/loop attachment on the side of the belt
- f) To turn off the belt, depress and hold the Start/Stop key(6) for 2 seconds.

If need be, the length of the belt can be extended by using the hook/loop attachments on the belt extender (4).

To turn off the belt, depress and hold the Start/Stop key(6) for 2 seconds and the LED which indicated the level of warmth will turn off.

5 Safety precautions

- Read all instructions carefully
- Do not puncture the warming belt with any pins, needles or other pointed objects.
- Warming Belt
 - do not use while folded or wrinkled
 - do not intended for use by sensation impaired people or by very small children and animals
- Use by minors is recommended only under direct supervision
Individuals who have impaired perception of heat, should only use the belt under supervision.
- Do not fold the belt while still warm
- Do not lay heavy or pointed objects on the belt so as not to damage the internal wires and electronics
- Do not expose the battery/controller to water or moisture
- Do not expose the belt to other heat producing objects during use.

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6 Product contains

- 1) Rechargeable warming belt
- 2) Controlling electronics(controller) and Lithium Ion rechargeable battery
- 3) Pocket for controller with hook and loop closure
- 4) Belt extender

7 Technical data

Model/Type	LH-115
Charger	100-240V~ 50-60Hz/8,4V
Product dimensions	ca. 100cm x 16,5cm
Weight	ca. 250g

For the sake of the environment, we request that you take packaging material to be recycled (collection points).

Disposal of appliances and batteries

After the implementation of the European Directive 2002/96/EU in the national legal

system, the following applies:

Electrical and electronic devices may not be disposed of with domestic waste.

Consumers are obliged by law to return electrical and electronic devices at the end of their service lives to the public collecting points set up for this purpose or point of sale. Details to this are defined by the national law of the

respective country. This symbol on the product, the instruction manual or the package indicates that a product is subject to these regulations. By recycling, reusing the materials or other forms of utilising old devices, you are making an important contribution to protecting our environment.

8 Guarantee

The guarantee period is 12 months from the date of purchase and covers all manufacturers' errors in material and quality. This guarantee only applies if the instructions for use have been followed and becomes ineffective if force has been applied to the appliance or it has been used in any but the correct and proper manner, or if it has been serviced by any unauthorised person.

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