**GB Instruction Manual**

**Item no. LH-885B**

**Relax Wrap**

**Please read before using!**

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**Item no. LH-885B**

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# English

## Introduction




Thank you for choosing the Relax Wrap. The Relax Wrap is a product of high standard that will last a long time if properly used and handled. There are four heat settings level, the Relax Wrap has a safety switch that will switch off automatically after 120 minutes.

The heat is regulated to ensure a permanent even heat distribution. If in use you will not experience the changes between on and off, or heating up and cooling down. A sensor cable provides for double protection against overheating.

The product is tested according to European standard and therefore technically safe. Please read the safety advice in your user manual.

## Relax Wrap Components

1. Relax Wrap
2. High Heat with red LED indicator
3. Medium Heat with orange LED indicator
4. Low Heat with yellow LED indicator
5. Warm Heat with green LED indicator
6. ON/OFF button
7. Connection clutch
8. Controller

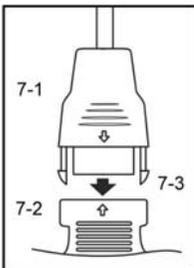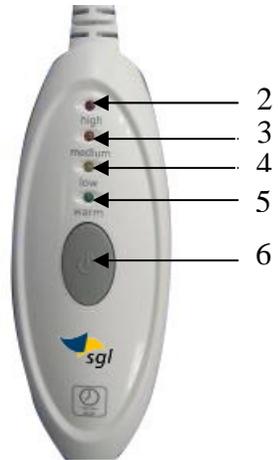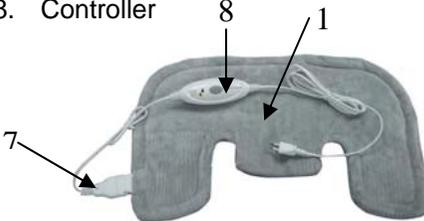- 7-1. Controller's Connector  
7-2. Wrap's Connector  
7-3. Locking Clips

**The appliance is not to be used by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction. Children being supervised not to play with the appliance.**

## **HOW TO USE**

Relax Wrap is specifically designed to fit the neck and shoulder area, where users experience the most physical signs of tension and stress. Drape Relax Wrap comfortably on your shoulders making sure that the neck flap is properly in place. Weighted edges and inter-grated magnets help to secure Relax Wrap in place. Note: This Relax Wrap utilizes magnets to fasten at the neck area. Remove necklaces and other metal accessories that could become attached to these magnets.

Make sure Connection Clutch (7) is securely connected.

The Relax Wrap is used for well-being and not for any medical purposes. Heat influences the blood flow and relaxes the muscles. The Wrap itself can be removed for washing.

### **Turn on heat**

Put plug into 220V AC electrical outlet. The Relax Wrap is ready for use. Select from 4 heat settings by pressing ON/OFF button successively to achieve desired heat level. To turn off, you must scroll through the remaining heat settings by pressing the ON/OFF button sequentially. The 4 heat settings with different LED indicators:

- The green LED indicated for Warm level.
- The yellow LED indicated for Lower level.
- The orange LED indicated for Medium level.
- The red LED indicated for High level.

Note: If the heat level indicators do not light up or stay on, the plug or control connector may not be firmly and completely connected. Unplug from wall outlet, check the Wrap connection again and re-plug into outlet.

### **Automatic switch off**

120 minutes after switching on, the Relax Wrap will switch off automatically. If the Relax Wrap should stay in use, you have to press ON/OFF button again.

**Attention:** Extended use may lead to burned skin.

## Using Moist Heat with Relax Wrap:


In some cases, you may find that moist heat provides more comfort and tension relief. When using moist heat, you should protect furniture and clothing from becoming damp.

1. Dampen one side of the Relax Wrap with water by patting the surface with a wet sponge or cloth or use a spray bottle to apply a fine mist. Do not thoroughly soak Wrap.
2. Do not spray or dampen the connector as per drawing 2.

**DANGER:** Do not immerse Wrap or connector in water or place under running water.


Failure to comply may result in product failure and/or personal injury.

Note: The electrical connector must be thoroughly dry before establishing electrical connection.
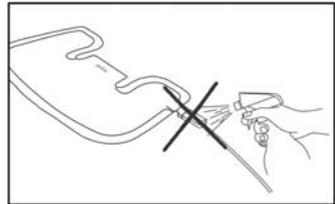
## Maintenance and Storage

The Relax Wrap can be washed in the washing machine after removing the Controller. Remove the mains cable at the connection clutch (7). You may return the Relax Wrap either to its shipping carton or store flat in a cool, dry place. Take care not to damage the surface of the Relax Wrap with sharp objects. The control unit must not be allowed to wash.
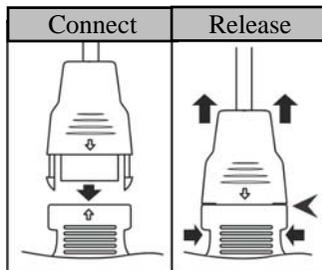
# Symbol explanation

Please, read the instructions

Do not use if folded or creased

Do not insert pins or needles

## Important Safety hints RETAIN MANUAL FOR FURTHER USE

- Read the user manual thoroughly
- Retain the manual for future use
- Do not prick in with needles or other sharp objects.
- Relax Wrap
  - only plug into adequate power point
  - do not use if wet or damp
  - do not use on heat insensitive persons, children or pets
  - do not use if bent or folded
  - do not use on helpless persons
  - do not fall asleep if Relax Wrap is switched on
  - do not fold Relax Wrap if still warm
- do not place heavy or edged objects on the Relax Wrap, which may damage or bend the electrical wires inside the Relax Wrap (eg. Suit cases )

- Keep plug and connecting cable dry at all times.
- The switch will get warm and therefore should not be covered, danger of overheating.
- Persons who only experience heat to a certain level should only use the Relax Wrap under supervision
- If the Relax Wrap or parts of it are damaged the Relax Wrap should not be used
- After use pull plug out off power point
- Do not place other heaters or heat producing objects onto the switched on cushion.
- Always supervise if used for children
- Avoid creasing the appliance
- Examine the appliance frequently for signs of wear or damage. If there are such signs or if the appliance has been misused, return it to the supplier prior to any further use
- This appliance is not intended for medical use in hospitals
- If the supply cords damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- Do not use if wet.
- When storing the appliance, allow it to cool down before folding.
- Do not crease the appliance by placing items on top of it during storage.
- When the heating cushion is likely to be used for an extended period, for instance by the user falling asleep, controls should be adjusted to a setting recommended for continuous use.
- The detachable control unit is only to be used with the types that are marked on the appliance.

## Delivery contents

- Relax Wrap
- User manual

## Technical data

Model:	LH-885B
Alternating current	220-240V~ 50Hz
Max. power	ca. 100W
Measurements	550*330mm

**For the sake of the environment**, we request that you take packaging material to be recycled (collection points).

### **Disposal of appliances and batteries**

After the implementation of the European Directive 2002/96/EU in the national legal system, the following applies:

Electrical and electronic devices may not be disposed of with domestic waste. Consumers are obliged by law to return electrical and electronic devices at the end of their service lives to the public collecting points set up for this purpose or point of sale . Details to this are defined by the national law of the respective country. This symbol on the product, the instruction manual or the package indicates that a product is subject to these regulations. By recycling, reusing the materials or other forms of utilising old devices, you are making an important contribution to protecting our environment.

### **Guarantee**

The guarantee period is 12 months from the date of purchase and covers all manufacturers' errors in material and quality. This guarantee only applies if the instructions for use have been followed and becomes ineffective if force has been applied to the appliance or it has been used in any but the correct and proper manner, or if it has been serviced by any unauthorised person.

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