




無線滑動熱力按摩墊
Gliding Massage Heated Cushion

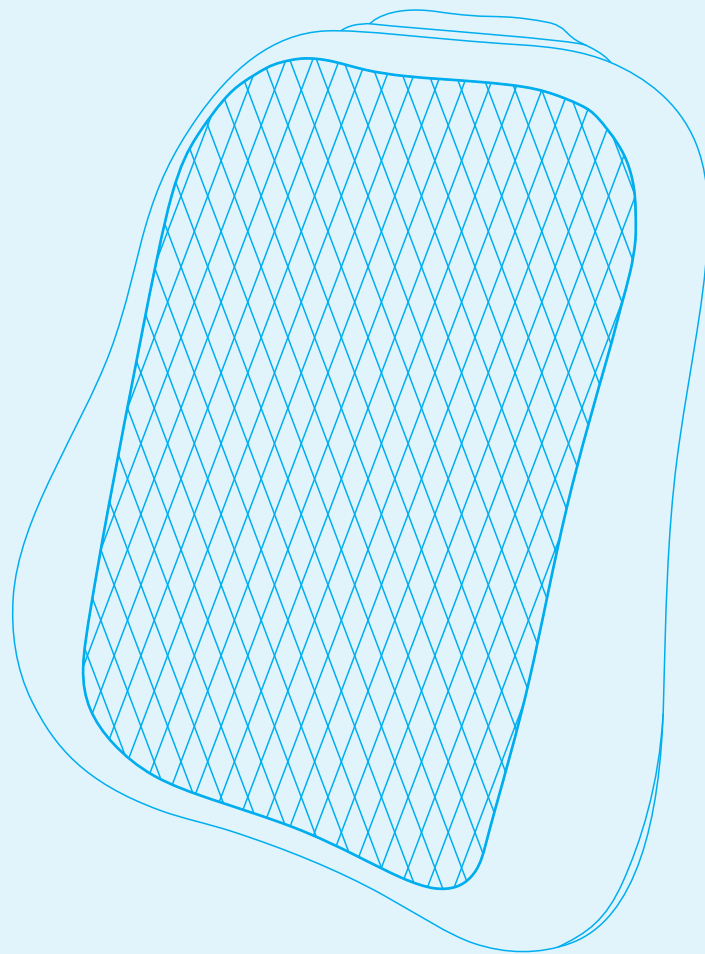

MS812
使用說明書
Instruction Manual

- 使用產品前，請先細閱本使用說明書，並請妥善保存以便將來翻閱。
- 本說明書及其修正版本（如有）可於本公司網站www.origin.hk下載。
- 本說明書所載的功能、外觀、圖片及數據皆以實物為準。
- Read all instructions before using the product. Keep this manual for future reference.
- This instruction manual is also available at www.origin.hk.
- In case of any discrepancies between this manual and the actual product—such as functions, drawings, photos, or data—the actual product shall prevail.

- 使用按摩墊前，請先細閱使用說明書。
Read all instructions carefully before using the Cushion.
- 必須按照使用說明書指示操作本按摩墊，否則可能引致產品損壞、短路、火警或人命損傷。
Use this Cushion only as described in this manual. Any other use not recommended by the manufacturer may cause fire, electrical shock or injury to persons.
- 打開包裝，檢查按摩墊。如發現產品有損壞，請勿使用及即時聯絡供應商。
Open the package and ensure that the Cushion is undamaged. In case of doubt, do not use the Cushion and contact your supplier.
- 包裝物料並非玩具，不可讓小孩玩耍。
Do not let children play parts of package.
- 為防止損害充電線，請勿把充電線掛放在檯邊上或接觸發熱的物體。
Do not hang the charging cable on the edge of table or touch the surface with high temperature in order to avoid damage to the power cord.
- 使用時，切勿把本產品放置於火爐或發熱的物體附近。
This product cannot be used beside the stove and other high temperature objects.
- 本機不是為體質、感覺或智力有障礙人士(包括小孩)或缺乏經驗及知識者使用，除非他們曾經接受負責安全人士的監督或指示，大人應指導小孩不可把玩電器。
This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- 本產品內有鋰電池，禁止拆解、撞擊、擠壓或投入火中。若出現嚴重鼓脹，請勿繼續使用。請勿置於高溫環境。
This product contains a lithium battery, and it is for hidden to disassemble, hit, squeeze or throw it into fire. Do not use if severe swelling occurs. Do not place in a high temperature environment.
- 每存放三個月充電一次，避免電池損壞。
Make a charging when not use the product for three months to avoid damage of battery.
- 在充電期間，切勿離開本產品。
Do not leave the appliance unattended while it is charging.

- 本產品內置充電電池，請使用本產品提供的USB充電線替產品進行完全充電。
This product is driven by batteries. Please use the USB charging cable provided to fully charge the product before use.
- 充滿電後，必須拔除充電線。長時間充電，可能造成充電電池過熱及燃燒等危險。
Make sure to unplug the charging cable after fully charged. Charging for a long time may cause hazards such as rechargeable battery overheat and burn.
- 不可在潮濕或有積水的環境使用，不可將本機放於浴缸或大型容器上面。
The appliance is not intended for use in wet or damp locations. Never locate the appliance where it may fall into a bathtub or other water container.
- 切勿於高溫表面上或附近使用本產品。
DO NOT OPERATE this appliance on or near any hot surface.
- 本產品並非玩具，請把本產品放置在兒童不能觸及的地方，以免孩童拉扯充電線或推翻本產品，發生意外。
Extreme caution is necessary when the Cushion is used by or near children or invalids and whenever the Cushion is left operating unattended.
- 本按摩墊不適用於室外及作工業用途。
This Cushion is intended for household use ONLY and not for industrial or outdoor use.
- 請勿在出氣口阻塞的狀況下操作本按摩墊。請保持出氣口通暢，避免被棉屑，毛髮等物阻塞。
Do not operate the Cushion with the air openings blocked. Keep the air openings free of lint, hair, and the likes.
- 不慎將液體傾倒在按摩墊上時，請立即將充電線從電源插座中拔除，以免導致觸電。
Remove charging cable from the electrical outlet immediately should liquid be accidentally spilled onto the Cushion as this may result in electric shock.
- 請勿在使用噴霧產品時或於氧氣不足的場所操作本按摩墊。
Do not operate where aerosol (spray) Cushion are being used or where oxygen is being administered.
- 請勿將任何物體或易燃材料掉落或插入按摩墊的任何開口中。這可能會導致觸電、故障或火災風險。
Do not drop or insert any object or flammable material into any openings of the Cushion. This may result in electric shock, malfunction or risk of fire.

重要安全提示

Important Safety Information

無線滑動熱力按摩墊
Gliding Massage Heated Cushion **BRIST**

- 請勿在室溫超過攝氏40度以上使用按摩墊。
Do not use the Cushion in temperatures higher than 40°C.
- 操作過程中不要蓋住按摩墊。這可能導致摩打溫度升高以及隨後的故障或火災。
Do not cover the Cushion during operation. This may result in a rise in motor temperature and subsequent malfunction or fire.
- 閃電情況發生時請勿使用按摩墊。
Do not use the Cushion in the event of a lightning storm.
- 若是將按摩墊由低溫的儲存處移至溫暖環境使用時，可能會因為機件內水氣聚集而影響本機正常運作，甚或造成機件故障。建議移動後，先將按摩墊在正常室溫放置1小時後再開始使用。
If the Cushion is brought into a warm environment from a cold place where it was stored, wait for an hour before use as there may be water droplets on mechanical parts due to condensation. Using the Cushion under such condition may result in malfunction.
- 不可站在按摩墊上，必須坐在座位上才能使用本按摩墊。
Do not stand on or in Cushion. Use only while seated.
- 使用按摩墊時請勿入睡。
Do not fall asleep when using the Cushion.
- 請避免在按摩前/使用產品前的1小時內進食。
Avoid food intake for at least an hour before massage / usage.
- 若產品的布料織物撕裂且內部零件暴露，請勿使用本按摩墊。
Do not use the Cushion if the fabric of the Cushion is torn and the internal components are exposed.
- 請勿拆卸、掉落或使產品受到衝擊或振動。
Do not dismantle, drop, or subject the product to shock or vibration.
- 為避免按摩效果減低或造成傷害，請勿將按摩墊與其它具有療效之設備或電暖墊等同時使用。
Do not use the Cushion in combination with other therapeutic equipment or electric blanket, etc, as this may result in ineffectiveness or injury.
- 請勿將重物堆放於按摩墊，以免造成意外傷害或機件故障。
Do not rest or place heavy objects on the Cushion as this may result in malfunction.
- 在使用中，可能會聽到按摩墊發出一些噪音。這是由於產品的結構，這是正常現象。
Some noise may be heard coming from the Cushion while in use. This is due to the structure of the product and is normal.
- 請勿在飲酒後或身體感到任何不適時使用按摩墊。
Please do not use the Cushion if you are not feeling well or under the influence of alcohol.
- 產品在啟動時會產生熱度。對熱度不敏感的人使用產品時必須小心謹慎。
The product has a heated surface. User who are insensitive to heat must be careful when using the product.

- 請勿擅自拆卸、修理或改造產品。碰觸內部機件可能導致機件故障或觸電意外。
Do not disassemble the cover of the product. Touching the internal components of the product may result in malfunction or an electric shock.
- 按摩時若身體感到任何不適，請立即停止操作。身體產生不正常或嚴重疼痛時，請尋求醫療諮詢。
Stop using immediately should you experience discomfort and consult your doctor if you experience any abnormal or severe pain.
- 為了達到最佳效果，每個自動按摩程式的特設時長為10分鐘。長時間連續使用可能會觸發過熱保護裝置，因此再次使用前請確保產品休息10分鐘。
Each auto massage program is designed to last 10 minutes for optimal effectiveness. Ensure that the product rests for 10 minutes before it is being used again as continuous usage over a long period may activate the overheat protection device.
- 建議每日使用時間請勿超過約30分鐘。過多的按摩可能會造成刺激過度而導致反效果。
It is recommended that your usage do not exceed a total of approximately 30 minutes per day as this may result in excessive stimulation and adverse effects.
- 請勿在毯子或枕頭下操作本產品。過熱可能造成火災、觸電或人員受傷。
Do not operate the product under the blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 如果您正在接受任何藥物治療或屬於以下任何類別，請在使用前尋求醫療諮詢：
心臟病、高血壓或低血壓、深部靜脈血、懷孕、糖尿病、體內植入了醫療器材(例如心臟節律器、自動心臟復律去顫器)、最近動過手術、患有任何骨骼疾病，如骨折、骨質疏鬆。
If you are under any medication or fall into any of the following categories, please consult your health care provider first before using:
Heart diseases, High, or low blood pressure, Deep Vein Thrombosis (DVT), Pregnant, Diabetes, Individuals fitted with medical implants (e.g. pacemaker, automatic cardiac defibrillator), Recently underwent surgical operations, Those with any bone conditions like fractures, osteoporosis.
- 本產品非自我診斷或治療之用；若需診斷或治療時，應由專業醫療人員處理。
This product is not intended for self-diagnosis or self-treatment of conditions that should be managed by a qualified health care provider.
- 定期按摩有助於身心健康。請適度的按摩，以避免過度刺激。保持均衡的按摩習慣，才能充分體驗按摩帶來的益處。
For best results, enjoy the massage regularly, while keeping in mind that moderation is key to avoid excessive stimulation and ensure a balanced experience.

- 必須使用本產品附配的USB線為本機充電。
Use only the supplied USB cable to charge the product.
- 依規格表列使用相應的充電器。
Use the corresponding charger according to the specification table.
- 為避免過熱和損壞裝置或電池，切勿在充電時使用本產品。
Do not charge while operating the product.
- 充滿電後，必須拔除充電線。長時間充電，可能造成充電池過熱及燃燒等危險。
To avoid overheat of rechargeable batteries, remove the USB cable after fully charged.
- 切勿在USB充電線纏結的情況下為產品充電。
Do not charge the product when the USB cable is tangle.
- 充電前，須確保產品的電源已關閉。
When charging, make sure that the power is off.
- 切勿使用有殘缺或鬆動的USB充電器。使用前確保USB充電器及充電線已充分抹乾及使用時插穩。
Never use deficient or loosen USB adapter. Make sure the USB charger and charging cable is dry before use and inserted firmly when use.
- 切勿使用無線充電。
Never use wireless charging.
- 充電期間，切勿離開本產品。
Do not leave the appliance unattended while it is charging.
- 充電進行中，充電器及熱力按摩墊會變熱，此乃正常現象。(充電器不包括在內)
The adaptor and Massage Heated Cushion will become hot while charging, which is normal. (Adaptor not included)
- 切勿在USB充電器、充電線及產品上放重物或讓其接觸灼熱的表面；切勿夾住、改裝、用力拉扯或扭曲USB充電器、充電線及產品。
Do not put heavy objects on the charger cable, USB adaptor and appliance or let it touch any hot surface. Never clip, modify, stretch or twist the USB charger, USB charging cable and appliance.
- 切勿在連接產品的充電線上施加任何壓力，充電線可能因此造成磨損和斷裂。
Do not put any stress on the charging cable where it connects to the appliance, as the cable may fray and break.
- 注意：持續使用本產品會使電池溫度升高，此乃正常現象。須待產品冷卻後再充電。
Continuously using the appliance will raise the temperature of the battery pack, which is normal. Wait until the product cool down before recharging.
- 請勿用金屬等接觸連接器而使其短路。
Do not short-circuit by allowing metals or other materials to come into contact with the socket terminals.


- 如本產品長期不使用時，建議每三個月充電一次，避免電池損壞導致產品不能使用。
It is recommended to charge it every 3 months to avoid battery damage and affect product functions if the product is not working for a long time.
- 充電及使用時間視乎設定和環境因素而定；實際結果可能有所差異。
Charging and operation time vary with settings and environmental factors; actual results will vary.
- 本產品包含只能由技術人員更換的電池。客戶不可自行更換。
The appliance contains batteries that are only replaced by skilled persons. The battery which is not self-replaceable by the user.
- 當發現電池漏液時，應立即停止使用產品及切勿充電，並聯絡服務中心進行檢查及維修。
When find the leaking batteries, stops use/charge the appliance immediately. Contact service centre for assistance.
- 若本產品出現嚴重膨脹，請立即停止使用，並聯絡服務中心安排檢查及維修。
If the appliance has severe swelling, please stop using it immediately and contact service centre for checking and repairing.
- 切勿毀壞或焚燒電池組，因為它們在高溫下會爆炸。它們必須根據當地法律進行處置。
Do not mutilate or incinerate the battery packs, as they will explode at high temperatures. They must be disposed of according to local laws.
- 為避免影響電池壽命，使用時因電量耗盡而自動停止時，請勿嘗試再次啟動，也請勿放置不管，以免電池過度放電。
To avoid affecting the battery life, when the product stops automatically due to the exhausted battery during use, do not attempt to restart it or leave it alone to avoid over discharge.
- 切勿將產品存放在室外或任何車輛中，因過熱會導致電池爆炸。
Do not store the appliance anywhere outdoors or in any vehicle as excessive heat can cause the batteries to explode.
- 本產品中的充電電池為消耗品。當即使正確充電後，使用時間仍明顯變短，電池壽命可能到了。
The rechargeable battery in this product is a consumable product. When the service time is still significantly shorter even after correct charging, the battery life may reach the end of its service life.
- 本產品廢棄時請取出電池，並交給官方指定回收站，切勿與一般生活垃圾堆放一起，以免污染環境。
Take out the battery before disposal and take it to the official designated recycle bin. Do not stack it with general domestic waste to avoid polluting the environment.
- 產品長期不使用時，請在常溫低濕度的環境下儲存，以免不良的環境導致產品功能損壞。
When the product is not working for a long time, please make sure it can storage under low humidity or room temperature in case of environment factors to the product function damage.

- 在大約5°C至35°C的溫度下充電和儲存本產品。
Charge and store this appliance at approximately 5°C to 35°C.
- 它被分類為第 III 類電器。III 類電器必須僅使用安全特低電壓 (SELV) 供電，除非它們是電池供電的，且電池是一次性電池或在電器外部充電的二次電池。除非它們是電池供電的，並且電池是一次性電池或在設備外部充電的二次電池。
It is classified as a class III appliance. Class III appliances must be supplied exclusively with Safety Extra Low Voltage (SELV), unless they are battery-operated and the battery is either a primary battery or a secondary battery charged outside the appliance.
- 切勿在陽光直射、高溫(靠近火源、火爐、烈日等)和易燃氣體的地方充電或閒置產品。否則可能導致產品變形、變色、故障、爆炸或火災。
Do not charge or leave in the places exposed to direct sunlight, places with high temperatures (near fires, stoves, under the scorching sun, etc) and flammable gas. It may cause deformation, discoloration, malfunction, explosion or fire of the product.
- 為避免充電對電視機、電話等產生干擾(如產生雜音等)，充電時請距離電視機、電話等1米以上。
In order to avoid interference with TV sets, telephones, etc. (Such as noise), please keep more than 1 m away from these appliances when charging.
-  此產品需要一個特定的獨立電源單元將本產品連接到電源。A specific separated power supply unit is required to connect this appliance to power supply.
- 本產品可運行時間的長短會因電池溫度及其他使用條件而有所差異。
The operation time can run varies depending on the temperature of the battery and other usage conditions.
- 當室溫低於 5°C 或高於 35°C時，有可能無法充電。
The battery may not be charged when the room temperature is below 5°C or above 35°C.
- 請勿使用本產品配置的充電線為其他產品充電。
Do not use the USB charging cable configured with this appliance to charge other appliance.
- 注意：為了給電池充電，只能使用本產品隨附的充電線。
WARNING: For the purpose of recharging the battery, only use the USB charging cable that provided with this appliance.
-  本產品含有可回收材料，不可把本產品作為未分類廢物。請聯繫當地就近指定回收點。
 This product contains recyclable materials. Do not dispose this product as unsorted municipal waste. Please contact your local municipality for the nearest collection point.

各部名稱 Main Parts

無線滑動熱力按摩墊
Gliding Massage Heated Cushion

BRIST

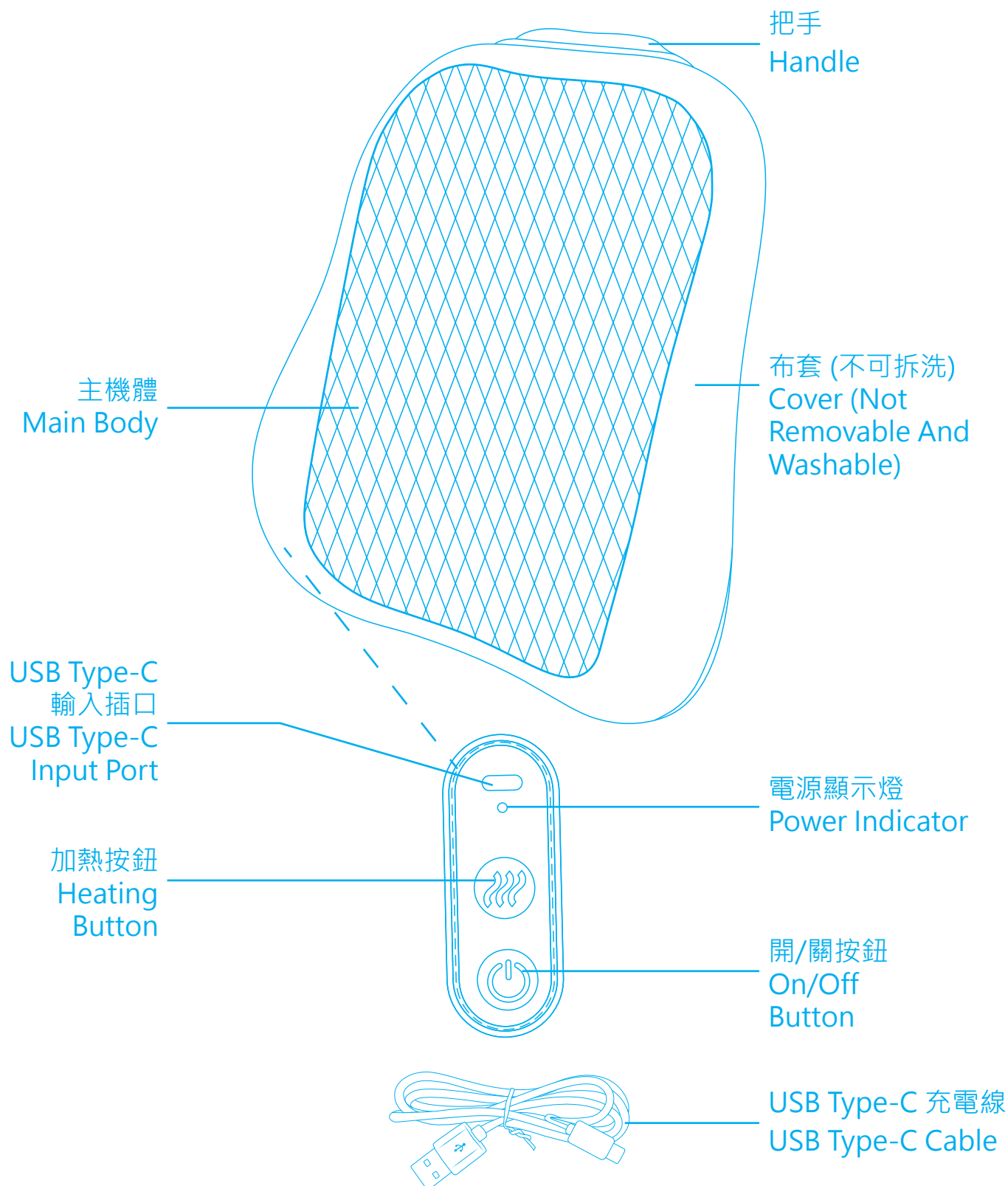

圖片只供參考，一切以實物為準
The picture is for reference only,
everything in kind shall prevail.

熱力按摩墊在日間坐在梳化時可放在背部使用，令你的背部變得溫暖，加上按摩功能正好促進血液循環及舒緩背部疲勞。本產品的設計不適合寵物使用。不正確使用會破壞熱力按摩墊的表面及發熱線的裝置。如有損毀，應立即停止使用，並交由指定維修服務中心檢查及修理。熱力按摩墊操作10分鐘後自動關掉。熱力按摩墊不同於發熱枕頭，另外手掌不及其他身體部份對熱的感應，用手觸碰是比較難測定熱力。

The Massage Heated Cushion can be placed on the back while sitting on the sofa during the day to warm your back. The warmth and massage function promote blood circulation and relieve back fatigue. This product is not designed for pets. Improper use can damage the surface of the Massage Heated Cushion and the heating element. If damaged, discontinue use immediately and take it to a designated service center for inspection and repair. The Massage Heated Cushion automatically turns off after 10 minutes of operation. Unlike a heating pillow, your palms are less sensitive to heat than other parts of the body, making it more difficult to sense the heat by touch.

按摩墊可使用背部按鈕操控。

The Cushion may be operated by the button on the back.

1. 長按開/關按鈕  啟動按摩墊，加熱功能同時開啟。

Long press the on/off button  to start Cushion. The heating function will be on.

2. 雙擊開/關按鈕  來改變按摩模式。(本機設有全背按摩和局部按摩兩種模式。)

Double press the on/off button  to change the massage modes. (There are two massage modes: full back massage mode and partial back massage mode.)

3. 按加熱按鈕 ，可開啟或關閉加熱功能。

Press the heating button  to turn the heating function on or off.

4. 再長按開/關按鈕 ，按摩墊關閉。

Long press the on/off button  again to turn off the Cushion.

一個按摩程序持續10分鐘後將自動關機。

The Cushion features an automatic shut-off function, which turns off automatically 10 minutes after being switched on.

自動關機

使用10分鐘後，熱力按摩墊會自動關機。如欲繼續使用熱力按摩墊，請重按開/關按鈕，熱力按摩墊將繼續另一次的10分鐘運作。

建議每次使用時間約10分鐘。確保再次使用按摩墊前休息10分鐘，因為長時間連續使用可能會啟動過熱保護設備。

注意: 過度連續長期使用，有可能灼傷皮膚。

Automatic Power Off

The Heated Cushion will automatically shut off after 10 minutes of use. To continue using the Heated Cushion, press the on/off button again. The Heated Cushion will continue to operate for another 10 minutes.

The suggested usage duration is 10 minutes. Let the Cushion rest for 10 minutes before it is being used again as continuous use over a long period may overheat device.

Caution: Excessive and prolonged use may cause burns.

使用完畢及待它冷卻後，請將熱力按摩墊放回包裝內。

After usage, put the Massage Heated Cushion back into the package.

保養

內部發熱元件及按摩組件不可放入水中清洗，只可使用濕軟毛巾抹擦表面污垢。切勿使用硬物刮花或破損按摩墊表面，以免發生危險。

清洗布套: 清洗前將USB Type-C充電線拔除，布套不可放入水中清洗，只可使用濕軟毛巾抹擦表面污垢。清洗後放於清涼地方涼乾，不可給陽光直接照射。

Maintenance

The internal heating element and massager assembly should not be washed in water. Use only a damp, soft towel to wipe any dirt from the surface. Avoid scratching or damaging the surface of the Cushion with hard objects to avoid danger.

Cleaning a cover: Before washing, unplug the USB type-C cable, the cover should not be washed in water. Use only a damp, soft towel to wipe any dirt from the surface. After washing, place it in a cool place to dry, away from direct sunlight.



不可乾洗及使用洗衣機清洗
Do not dry clean and machine wash



不可熨燙
Do not iron



不可使用機器烘乾
Do not tumble dry



請小心閱讀使用說明書
Read all the instructions



不可漂染
Do not bleach



不可用針或扣針戳穿按摩墊
Do not puncture the Cushion

1. 充電/使用時間：

充電時間：約4小時。

操作時間：最多約60分鐘(每節10分鐘)。

提示：根據電池壽命，充電電池的使用時間會持續縮短。

CHARGING/USING TIME:

Charging Time: Approx. 4 Hours.

Operation Time: Approx. 60 Minutes (10 mins/session).

Reminder : According to the battery nature, the using time of the rechargeable battery will be shortened continuously.

2. 充電: 將USB線連接到設備(例如:電源適配器(需自備))，另一端連接到按摩墊背面的電源輸入端口，為按摩墊充電。

CHARGING: Connect the USB cable to a device(e.g. power adaptor(not included)), and the power input port of the back of the Heated Cushion to charge the Heated Cushion.

- 首次充電時間: 最少充電約4小時。

First charging time: At least 4 Hours.

- 充電時充電指示燈會閃爍，充電完成後會長亮。

The charging indicator light will flash while charging and will stay on when charging is complete.

- 如本產品長期不使用時，建議每三個月充電一次，避免電池損壞導致產品不能使用。


It is recommended to charge it every 3 months to avoid battery damage and affect product functions if the product is not working for a longtime.

- 建議使用DC5V 2A電源充電，否則會使充電電池的壽命縮短。(電源適配器不包括在內)

It is recommended to use DC5V 2A power supply for charging to achieve the better battery performance. (Power adaptor not included)

- 充電及使用時間視乎設定和環境因素而定；實際結果可能有所差異。


Charging and operation time vary with settings and environmental factors; actual results will vary.

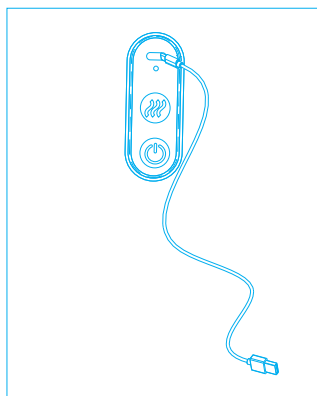

圖片只供參考，一切以實物為準
The picture is for reference only,
everything in kind shall prevail.

型號 Model Number	MS812
輸入 Input	DC5V 2A / USB-C Port
電池 Battery	Li-ion 2500mAh
充電 Charging	4小時/Hours(約/Approx.)
操作時間 Operation	10 分鐘mins/節session 60 分鐘/mins Max. (約/Approx.)
尺寸 Dimension	390(W) x 450(H) x 120(D)mm(約/Approx.)
淨重 Net Weight	1420克(gram) (約/Approx.)

源澤地圖 Origin Map

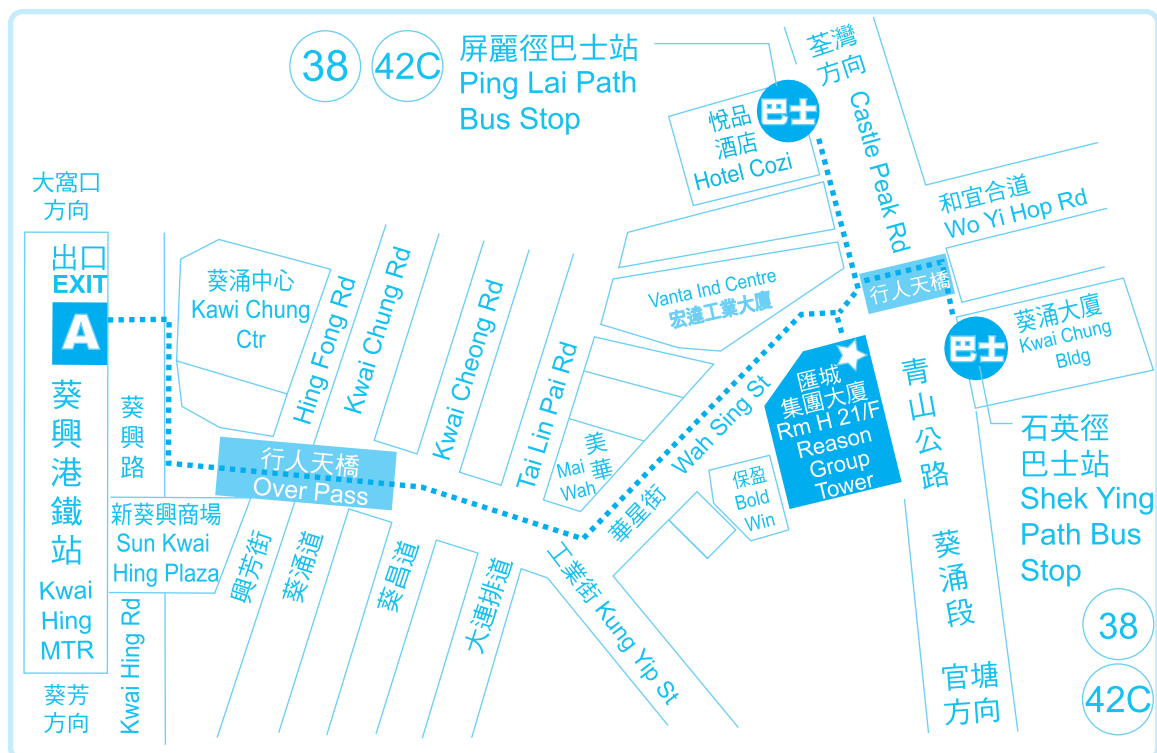

辦公時間

10:30am - 12:30pm (星期一至五)
02:30pm - 05:30pm (星期一至五)
10:30am - 12:30pm (星期六)
(星期日及公眾假期休息)

OFFICE HOURS

10:30am - 12:30pm (Mon to Fri)
02:30pm - 05:30pm (Mon to Fri)
10:30am - 12:30pm (Saturday)
(Closed on Sunday and Holiday)

保用

無線滑動熱力按摩墊
Gliding Massage Heated Cushion **BRIST**

產品：無線滑動熱力按摩墊

型號：MS812

此產品由購買日期起計算，本公司將根據下列條款及細則，提供一年保養服務(商業用途則為3個月)：

1. 於購買後一年內，經本公司技術員判斷為正常運作下損壞之產品，可獲免費維修和/或零件更換。

2. 以下情況不被包括在保養範圍內：

- 若產品是因為濫用、錯誤使用、疏忽及/或意外而遭致損壞；
- 未經本公司授權之自行處理、維修及/或改裝；
- 任何布套、USB Type-C 充電線，以及配件。

3. 需要檢查或維修時，請將產品送往本公司服務中心，並於修理妥當後自行取回。

4. 更換後之任何零件或配件，均屬本公司所有。

5. 本保養服務只適用於香港特別行政區。

6. 顧客必須同時出示零售商之正式發票，方可得到免費維修服務。

升級保用：(須於購買本產品後10天內完成以下步驟；升級保用不適用於商業用途)

A. 額外3個月保用：於本公司網站完成登記，可享共有15個月免費保用；或

B. 額外6個月保用：到訪並“追蹤”本公司 (Origin Marketing Limited) 的Facebook，即可享有共18個月免費保用。(如果已追蹤，則自動享有此升級保用，毋須再進行上述步驟)。

服務中心

源澤市務有限公司

香港新界葵涌青山公路403號(近華星街12號)匯城集團大廈21樓H室

www.origin.hk 電郵: info@origin.hk 電話: 2156 8238

www.origin.hk



Find us on



Instragarm



Warranty

Item : Gliding Massage Heated Cushion

Model : MS812

This product is warranted for 1 year (commercial use for 3 months) from the date of purchase under the following conditions:

1. For any defect, in the judgment of our technician, if it is caused under normal use, we are responsible for repairing the product, and/or replacing parts at free of charge.

2. This warranty DOES NOT apply to:

- If the product was damaged by abuse, misuse, negligence or accidents;
- Repair, alteration/modification by unauthorized technician;
- The cover, USB Type-C cable or accessories.

3. The customer is requested to bring the product to and from our Service Centre.

4. Any parts or accessories, which had been replaced, shall be of our property.

5. This warranty is valid in H.K.S.A.R. only.

6. Original invoice should be presented in order to get free service.

Extended Warranty: (Complete following procedure within 10 days from the date of purchase.

Extended warranty does not apply to commercial use)

A. Extra 3 months warranty: complete Online Warranty Registration on our website to enjoy a total of 15 months of warranty; or

B. Extra 6 months warranty: follow "Origin Marketing Limited" on Facebook to enjoy a total of 18 months of warranty. (If you have already followed, no need to do it again.)

Service Centre

ORIGIN MARKETING LIMITED

Unit H, 21/F., Reason Group Tower (nearby 12 Wah Sing Street),

403 Castle Peak Road, Kwai Chung, N.T., Hong Kong.

www.origin.hk Email: info@origin.hk Tel: 2156 8238

www.origin.hk



Find us on



Instragarm

