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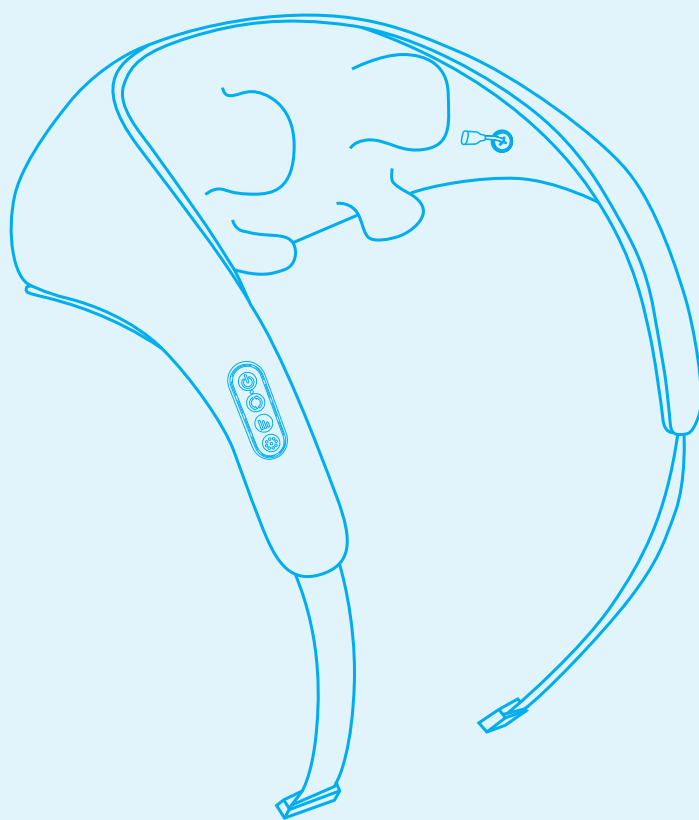
- 使用按摩器前，請先細閱使用說明書。
Read all instructions carefully before using the Massager.
- 必須按照使用說明書指示操作本按摩器，否則可能引致產品損壞、短路、火警或人命損傷。
Use this Massager only as described in this manual. Any other use not recommended by the manufacturer may cause fire, electrical shock or injury to persons.
- 打開包裝，檢查按摩器。如發現產品有損壞，請勿使用及即時聯絡供應商。
Open the package and ensure that the Massager is undamaged. In case of doubt, do not use the Massager and contact your supplier.
- 包裝物料並非玩具，不可讓小孩玩耍。
Do not let children play parts of package.
- 為防止損害充電線，請勿把充電線掛放在檯邊上或接觸發熱的物體。
Do not hang the charging cable on the edge of table or touch the surface with high temperature in order to avoid damage to the power cord.
- 使用時，切勿把本產品放置於火爐或發熱的物體附近。
This product cannot be used beside the stove and other high temperature objects.
- 本機不是為體質、感覺或智力有障礙人士(包括小孩)或缺乏經驗及知識者使用，除非他們曾經接受負責安全人士的監督或指示，大人應指導小孩不可把玩電器。
This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- 本產品內有鋰電池，禁止拆解、撞擊、擠壓或投入火中。若出現嚴重鼓脹，請勿繼續使用。請勿置於高溫環境。
This product contains a lithium battery, and it is for hidden to disassemble, hit, squeeze or throw it into fire. Do not use if severe swelling occurs. Do not place in a high temperature environment.
- 每存放三個月充電一次，避免電池損壞。
Make a charging when not use the product for three months to avoid damage of battery.
- 在充電期間，切勿離開本產品。
Do not leave the appliance unattended while it is charging.

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Important Safety Information

3D 指壓頸部熱力按摩器
3D Shiatsu Neck Heated Massager **BRIST**

- 本產品內置充電電池，請使用本產品提供的USB充電線替產品進行完全充電。
This product is driven by batteries. Please use the USB charging cable provided to fully charge the product before use.
- 充滿電後，必須拔除充電線。長時間充電，可能造成充電電池過熱及燃燒等危險。
Make sure to unplug the charging cable after fully charged. Charging for a long time may cause hazards such as rechargeable battery overheat and burn.
- 不可在潮濕或有積水的環境使用，不可將本機放於浴缸或大型容器上面。
The appliance is not intended for use in wet or damp locations. Never locate the appliance where it may fall into a bathtub or other water container.
- 切勿於高溫表面上或附近使用本產品。
DO NOT OPERATE this appliance on or near any hot surface.
- 本產品並非玩具，請把本產品放置在兒童不能觸及的地方，以免孩童拉扯充電線或推翻本產品，發生意外。
Extreme caution is necessary when the Massager is used by or near children or invalids and whenever the Massager is left operating unattended.
- 本按摩器不適用於室外及作工業用途。
This Massager is intended for household use ONLY and not for industrial or outdoor use.
- 請勿在出氣口阻塞的狀況下操作本按摩器。請保持出氣口通暢，避免被棉屑，毛髮等物阻塞。
Do not operate the Massager with the air openings blocked. Keep the air openings free of lint, hair, and the likes.
- 不慎將液體傾倒在按摩器上時，請立即將充電線從電源插座中拔除，以免導致觸電。
Remove charging cable from the electrical outlet immediately should liquid be accidentally spilled onto the Massager as this may result in electric shock.
- 請勿在使用噴霧產品時或於氧氣不足的場所操作本按摩器。
Do not operate where aerosol (spray) Massager are being used or where oxygen is being administered.
- 請勿將任何物體或易燃材料掉落或插入按摩器的任何開口中。這可能會導致觸電、故障或火災風險。
Do not drop or insert any object or flammable material into any openings of the Massager. This may result in electric shock, malfunction or risk of fire.

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Important Safety Information

3D 指壓頸部熱力按摩器
3D Shiatsu Neck Heated Massager **BRIST**

- 請勿擅自拆卸、修理或改造產品。碰觸內部機件可能導致機件故障或觸電意外。
Do not disassemble the cover of the product. Touching the internal components of the product may result in malfunction or an electric shock.
- 按摩時若身體感到任何不適，請立即停止操作。身體產生不正常或嚴重疼痛時，請尋求醫療諮詢。
Stop using immediately should you experience discomfort and consult your doctor if you experience any abnormal or severe pain.
- 為了達到最佳效果，每個自動按摩程式的特設時長為10分鐘。長時間連續使用可能會觸發過熱保護裝置，因此再次使用前請確保產品休息10分鐘。
Each auto massage program is designed to last 10 minutes for optimal effectiveness. Ensure that the product rests for 10 minutes before it is being used again as continuous usage over a long period may activate the overheat protection device.
- 建議每日使用時間請勿超過約30分鐘。過多的按摩可能會造成刺激過度而導致反效果。
It is recommended that your usage do not exceed a total of approximately 30 minutes per day as this may result in excessive stimulation and adverse effects.
- 請勿在毯子或枕頭下操作本產品。過熱可能造成火災、觸電或人員受傷。
Do not operate the product under the blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 如果您正在接受任何藥物治療或屬於以下任何類別，請在使用前尋求醫療諮詢：
心臟病、高血壓或低血壓、深部靜脈血、懷孕、糖尿病、體內植入了醫療器材(例如心臟節律器、自動心臟復律去顫器)、最近動過手術、患有任何骨骼疾病，如骨折、骨質疏鬆。
If you are under any medication or fall into any of the following categories, please consult your health care provider first before using:
Heart diseases, High, or low blood pressure, Deep Vein Thrombosis (DVT), Pregnant, Diabetes, Individuals fitted with medical implants (e.g. pacemaker, automatic cardiac defibrillator), Recently underwent surgical operations, Those with any bone conditions like fractures, osteoporosis.
- 本產品非自我診斷或治療之用；若需診斷或治療時，應由專業醫療人員處理。
This product is not intended for self-diagnosis or self-treatment of conditions that should be managed by a qualified health care provider.
- 定期按摩有助於身心健康。請適度的按摩，以避免過度刺激。保持均衡的按摩習慣，才能充分體驗按摩帶來的益處。
For best results, enjoy the massage regularly, while keeping in mind that moderation is key to avoid excessive stimulation and ensure a balanced experience.

- 必須使用本產品附配的USB線為本機充電。
Use only the supplied USB cable to charge the product.
- 依規格表列使用相應的充電器。
Use the corresponding charger according to the specification table.
- 為避免過熱和損壞裝置或電池，切勿在充電時使用本產品。
Do not charge while operating the product.
- 充滿電後，必須拔除充電線。長時間充電，可能造成充電池過熱及燃燒等危險。
To avoid overheat of rechargeable batteries, remove the USB cable after fully charged.
- 切勿在USB充電線纏結的情況下為產品充電。
Do not charge the product when the USB cable is tangle.
- 充電前，須確保產品的電源已關閉。
When charging, make sure that the power is off.
- 切勿使用有殘缺或鬆動的USB充電器。使用前確保USB充電器及充電線已充分抹乾及使用時插穩。
Never use deficient or loosen USB adapter. Make sure the USB charger and charging cable is dry before use and inserted firmly when use.
- 切勿使用無線充電。
Never use wireless charging.
- 充電期間，切勿離開本產品。
Do not leave the appliance unattended while it is charging.
- 充電進行中，充電器及熱力按摩器會變熱，此乃正常現象。(充電器不包括在內)
The adaptor and Heated Massager will become hot while charging, which is normal. (Adaptor not included)
- 切勿在USB充電器、充電線及產品上放重物或讓其接觸灼熱的表面；切勿夾住、改裝、用力拉扯或扭曲USB充電器、充電線及產品。
Do not put heavy objects on the charger cable, USB adaptor and appliance or let it touch any hot surface. Never clip, modify, stretch or twist the USB charger, USB charging cable and appliance.
- 切勿在連接產品的充電線上施加任何壓力，充電線可能因此造成磨損和斷裂。
Do not put any stress on the charging cable where it connects to the appliance, as the cable may fray and break.
- 注意：持續使用本產品會使電池溫度升高，此乃正常現象。須待產品冷卻後再充電。
Continuously using the appliance will raise the temperature of the battery pack, which is normal. Wait until the product cool down before recharging.
- 請勿用金屬等接觸連接器而使其短路。
Do not short-circuit by allowing metals or other materials to come into contact with the socket terminals.


- 如本產品長期不使用時，建議每三個月充電一次，避免電池損壞導致產品不能使用。
It is recommended to charge it every 3 months to avoid battery damage and affect product functions if the product is not working for a long time.
- 充電及使用時間視乎設定和環境因素而定；實際結果可能有所差異。
Charging and operation time vary with settings and environmental factors; actual results will vary.
- 本產品包含只能由技術人員更換的電池。客戶不可自行更換。
The appliance contains batteries that are only replaced by skilled persons. The battery which is not self-replaceable by the user.
- 當發現電池漏液時，應立即停止使用產品及切勿充電，並聯絡服務中心進行檢查及維修。
When find the leaking batteries, stops use/charge the appliance immediately. Contact service centre for assistance.
- 若本產品出現嚴重膨脹，請立即停止使用，並聯絡服務中心安排檢查及維修。
If the appliance has severe swelling, please stop using it immediately and contact service centre for checking and repairing.
- 切勿毀壞或焚燒電池組，因為它們在高溫下會爆炸。它們必須根據當地法律進行處置。
Do not mutilate or incinerate the battery packs, as they will explode at high temperatures. They must be disposed of according to local laws.
- 為避免影響電池壽命，使用時因電量耗盡而自動停止時，請勿嘗試再次啟動，也請勿放置不管，以免電池過度放電。
To avoid affecting the battery life, when the product stops automatically due to the exhausted battery during use, do not attempt to restart it or leave it alone to avoid over discharge.
- 切勿將產品存放在室外或任何車輛中，因過熱會導致電池爆炸。
Do not store the appliance anywhere outdoors or in any vehicle as excessive heat can cause the batteries to explode.
- 本產品中的充電電池為消耗品。當即使正確充電後，使用時間仍明顯變短，電池壽命可能到了。
The rechargeable battery in this product is a consumable product. When the service time is still significantly shorter even after correct charging, the battery life may reach the end of its service life.
- 本產品廢棄時請取出電池，並交給官方指定回收站，切勿與一般生活垃圾堆放一起，以免污染環境。
Take out the battery before disposal and take it to the official designated recycle bin. Do not stack it with general domestic waste to avoid polluting the environment.
- 產品長期不使用時，請在常溫低濕度的環境下儲存，以免不良的環境導致產品功能損壞。
When the product is not working for a long time, please make sure it can storage under low humidity or room temperature in case of environment factors to the product function damage.

- 在大約5°C至35°C的溫度下充電和儲存本產品。
Charge and store this appliance at approximately 5°C to 35°C.
- 它被分類為第 III 類電器。III 類電器必須僅使用安全特低電壓 (SELV) 供電，除非它們是電池供電的，且電池是一次性電池或在電器外部充電的二次電池。除非它們是電池供電的，並且電池是一次性電池或在設備外部充電的二次電池。
It is classified as a class III appliance. Class III appliances must be supplied exclusively with Safety Extra Low Voltage (SELV), unless they are battery-operated and the battery is either a primary battery or a secondary battery charged outside the appliance.
- 切勿在陽光直射、高溫(靠近火源、火爐、烈日等)和易燃氣體的地方充電或閒置產品。否則可能導致產品變形、變色、故障、爆炸或火災。
Do not charge or leave in the places exposed to direct sunlight, places with high temperatures (near fires, stoves, under the scorching sun, etc) and flammable gas. It may cause deformation, discoloration, malfunction, explosion or fire of the product.
- 為避免充電對電視機、電話等產生干擾(如產生雜音等)，充電時請距離電視機、電話等1米以上。
In order to avoid interference with TV sets, telephones, etc. (Such as noise), please keep more than 1 m away from these appliances when charging.
-  此產品需要一個特定的獨立電源單元將本產品連接到電源。A specific separated power supply unit is required to connect this appliance to power supply.
- 本產品可運行時間的長短會因電池溫度及其他使用條件而有所差異。
The operation time can run varies depending on the temperature of the battery and other usage conditions.
- 當室溫低於 5°C 或高於 35°C時，有可能無法充電。
The battery may not be charged when the room temperature is below 5°C or above 35°C.
- 請勿使用本產品配置的充電線為其他產品充電。
Do not use the USB charging cable configured with this appliance to charge other appliance.
- 注意：為了給電池充電，只能使用本產品隨附的充電線。
WARNING: For the purpose of recharging the battery, only use the USB charging cable that provided with this appliance.
-  本產品含有可回收材料，不可把本產品作為未分類廢物。請聯繫當地就近指定回收點。
 This product contains recyclable materials. Do not dispose this product as unsorted municipal waste. Please contact your local municipality for the nearest collection point.

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3D 指壓頸部熱力按摩器
3D Shiatsu Neck Heated Massager

BRIST

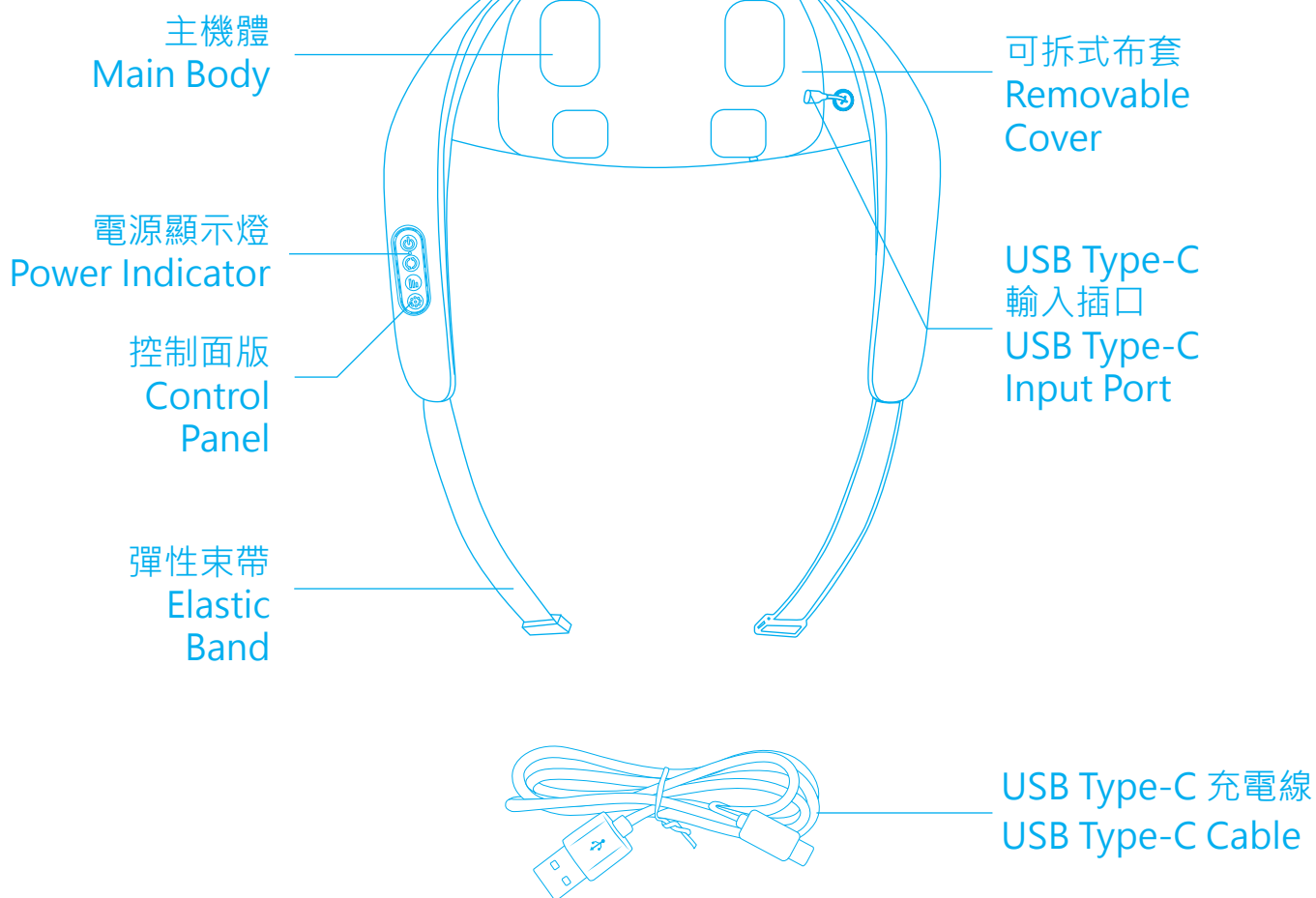圖片只供參考，一切以實物為準
The picture is for reference only,
everything in kind shall prevail.

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Notice: When using the product for the first time, start with low massaging speed for a short period to understand the features and determine your preferred intensity. Use the elastic tightening band to tie the device onto your body or hold the handles. Let the massage balls land evenly on your body.

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The Heated Massager can be placed on the neck or back or leg while sitting on the sofa during the day to warm your neck, back or leg. The warmth and massage function promote blood circulation and relieve neck, back or leg fatigue. This product is not designed for pets. Improper use can damage the surface of the Heated Massager and the heating element. If damaged, discontinue use immediately and take it to a designated service center for inspection and repair. The Heated Massager automatically turns off after 10 minutes of operation. Unlike a heating pillow, your palms are less sensitive to heat than other parts of the body, making it more difficult to sense the heat by touch.

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Choose the body area that you want to massage. The Massager is designed for neck, back or leg areas.

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Long press the on/off button  to start Massager.

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Press the  button to change the massaging direction.

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Press the  button to select two massaging speeds.

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Press the heating button , the Massager will start operation at Low Heat, press this button again set the Massager into High Heat or turn the heating function off.

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To turn off the Massager, press the on/off button  again.

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The Massager will turn off automatically 10 minutes after being switched on.

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Notice: Do not apply the massage directly on joints of elbows and knees, breasts, or abdomen.

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Automatic Power Off

The Heated Massager will automatically shut off after 10 minutes of use. To continue using the Heated Massager, press the on/off button again. The Heated Massager will continue to operate for another 10 minutes.

The suggested usage duration is 10 minutes. Let the Massager rest for 10 minutes before it is being used again as continuous use over a long period may overheat device.

Caution: Excessive and prolonged use may cause burns.

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After use, put the Heated Massager back into the package.

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Maintenance

The internal heating element and massager assembly should not be washed. Use only a damp, soft towel to wipe any dirt from the surface. Avoid scratching or damaging the surface of the Massager with hard objects as to avoid danger.

Cleaning a cover: Before washing, unplug the USB type-C cable, unzip the cover, and remove it. The cover can be washed with water. After washing, place it in a cool place to dry, away from direct sunlight.

Note: The cover is hand washable only, it must be completely dry before use.

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CHARGING/USING TIME:

Charging Time: Approx. 4 Hours.

Operation Time: Approx. 60 Minutes (10 mins/session).

Reminder : According to the battery nature, the using time of the rechargeable battery will be shortened continuously.

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CHARGING: Connect the USB cable to a device(e.g. power adaptor(not included)), and the power input port of the Heated Massager to charge the Heated Massager.

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The charging indicator light turn red while charging and the indicator light turns green when charging is complete.


- 如本產品長期不使用時，建議每三個月充電一次，避免電池損壞導致產品不能使用。
It is recommended to charge it every 3 months to avoid battery damage and affect product functions if the product is not working for a longtime.

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It is recommended to use DC5V 2A power supply for charging to achieve the better battery performance. (Power adaptor not included)

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
Charging and operation time vary with settings and environmental factors; actual results will vary.

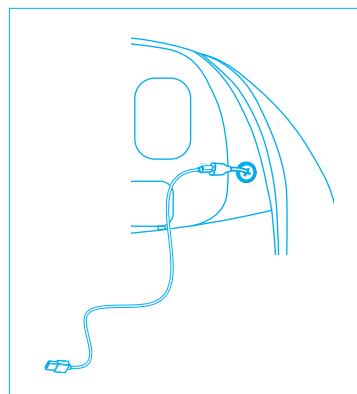圖片只供參考，一切以實物為準
The picture is for reference only,
everything in kind shall prevail.

型號 Model Number	MS813
輸入 Input	DC5V 2A / USB-C Port
電池 Battery	Li-ion 1800mAh
充電 Charging	4小時/Hours(約/Approx.)
操作時間 Operation	10 分鐘mins/節session 60 分鐘/mins Max. (約/Approx.)
尺寸 Dimension	260(W) x 160(H) x 140(D)mm(約/Approx.)
淨重 Net Weight	1170克(gram) (約/Approx.)

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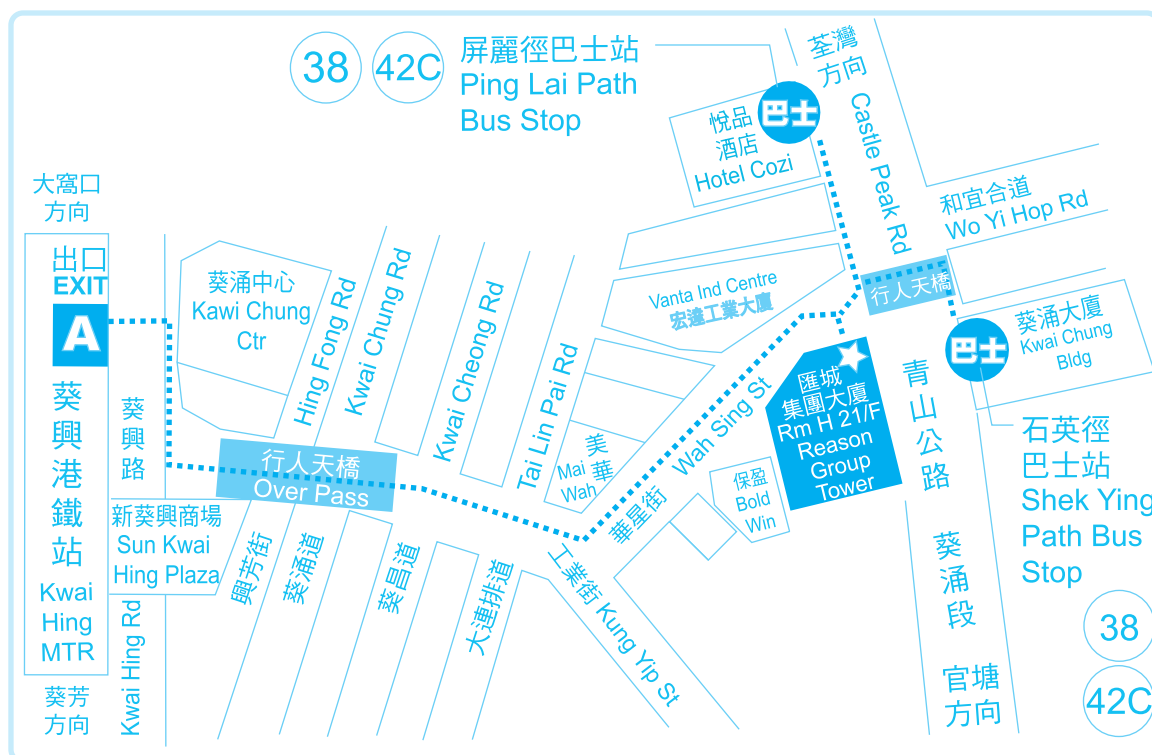OFFICE HOURS

10:30am - 12:30pm (Mon to Fri)
02:30pm - 05:30pm (Mon to Fri)
10:30am - 12:30pm (Saturday)
(Closed on Sunday and Holiday)

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Item : 3D Shiatsu Neck Heated Massager

Model : MS813

This product is warranted for 1 year (commercial use for 3 months) from the date of purchase under the following conditions:

1. For any defect, in the judgment of our technician, if it is caused under normal use, we are responsible for repairing the product, and/or replacing parts at free of charge.

2. This warranty DOES NOT apply to:

- If the product was damaged by abuse, misuse, negligence or accidents;
- Repair, alteration/modification by unauthorized technician;
- The cover, USB Type-C cable or accessories.

3. The customer is requested to bring the product to and from our Service Centre.

4. Any parts or accessories, which had been replaced, shall be of our property.

5. This warranty is valid in H.K.S.A.R. only.

6. Original invoice should be presented in order to get free service.

Extended Warranty: (Complete following procedure within 10 days from the date of purchase.

Extended warranty does not apply to commercial use)

A. Extra 3 months warranty: complete Online Warranty Registration on our website to enjoy a total of 15 months of warranty; or

B. Extra 6 months warranty: follow "Origin Marketing Limited" on Facebook to enjoy a total of 18 months of warranty. (If you have already followed, no need to do it again.)

Service Centre

ORIGIN MARKETING LIMITED

Unit H, 21/F., Reason Group Tower (nearby 12 Wah Sing Street),

403 Castle Peak Road, Kwai Chung, N.T., Hong Kong.

www.origin.hk Email: info@origin.hk Tel: 2156 8238

www.origin.hk

Find us on f

Instragarm